

Link-A-Light R-Phycoerythrin (R-PE) Conjugation Kit (3 labeling reactions)

Catalog No.:	AC004-003
Quantity:	3 x 100 µg
Reaction:	3 labeling reactions
Kit Contents:	AC004-003-M1: Link-A-Light Mix (Ready to Use) AC004-003-R1: Link-A-Light Modifier reagent (Ready to Use) AC004-003-Q1: Link-A-Light Quencher reagent (Ready to Use)
Conjugate:	PE
Background:	Link-A-Light is a one step antibody labeling kit, requiring just 30 seconds hands-on time. The antibody is covalently bonded to the label in a directional and controlled process at near-neutral pH.

The Link-A-Light conjugation kit allows PE conjugations to set up in seconds, simply by adding a solution of the antibody to be labeled to a proprietary lyophilised mixture containing PE (Figure 1).

By circumventing the desalting or dialysis steps that commonly interrupt traditional protein conjugation procedures, Link-A-Light technology can be used to label small quantities of protein for FACS analysis (see Principle).

The Link-A-Light conjugation kit allows PE conjugations to set up in seconds, simply by adding a solution of the antibody to be labeled to a proprietary lyophilised mixture containing PE (Figure 1).

Label: R-Phycoerythrin (R-PE) is a fluorescent protein from the phycobiliprotein family, and is isolated from red algae. The absorbance spectrum of R-PE covers a broad range of excitation wavelengths, which provides an advantage for multi-color immunofluorescent staining or cell sorting. R-PE is one of the most intensely fluorescent phycobiliprotein having orange fluorescence. It is significantly brighter and more photostable than conventional organic fluorophores and it has a high quantum yield.

Applications: Labeling

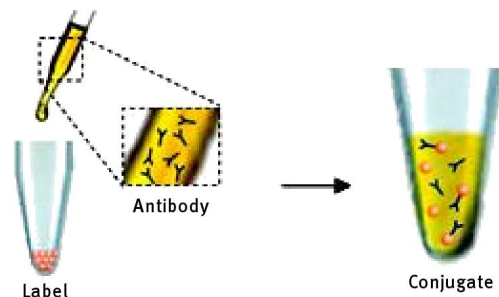


Figure 1

Principle of the Procedure:

Upon dissolution of Link-A-Light mixture with a solution of the antibody (or other biomolecule to be labeled) proprietary chemicals in the mixture become activated. This results in the *directional, covalent* bonding of the antibody to the Fluorescent label in a

gentle and controlled process at near-neutral pH. The hands-on time for the entire procedure is usually 20-30 seconds. Link-A-Light makes it possible to label antibodies and other biomolecules with PE with ease, and eliminates the need for secondary detection reagents in FACS experiments. Direct labeling of antibodies simplifies and improves data quality in multicolor experiments by eliminating problems caused by dissociation and crossover of secondary reagents.

Frequently asked Questions:

Q1. What functional groups do I need on my protein? Link-A-Light requires amine groups on the molecule to be labeled. Most proteins have lysine and/or alpha-amino groups. All antibodies will have multiple amine functions. **Q2.** Do I need to purify the conjugate? No. The chemicals used in Link-A-Light are deactivated by the Quencher, and the by-products are benign and do not need to be removed.

Materials Required but Not Provided:

Antibody to be Labeled: In view of the large size of PE (240kDa), the amount of antibody used in a labeling reaction must always be less than the pack size of Link-A-Light-PE, in order that the PE is in a slight molar excess. The best ratio for any new antibody reagent must be determined by experimentation but 50-60 µg of IgG antibody for every 100 µg of Link-A-Light PE usually gives optimal results. The 60 µg quantity corresponds to an Ab:PE molar ratio of 1:1. The volume in which the antibody is added ideally should be around 40 µl (100 µg pack size), and around 400 µl (1 mg pack size). Where the concentration of antibody is relatively low, and where it is impractical to concentrate the antibody, up to twice the volume stated above (i.e. 80 µl for the 100 µg PE pack size) may be added without any significant loss in conjugation efficiency. **Sample Buffer:** Ideally, the antibody to be labeled should be in 10-50mM amine-free buffer pH range 6.5 to 8.5. However, many buffers outside these limits of concentration and pH can be accommodated. Modest concentrations of Tris buffer are also tolerated. **Compatibility of Buffers and Buffer Additives:** Amine-free buffers, including MES, MOPS, HEPES and phosphate are compatible if they are in the concentration range 10-50mM and have pH values in the range 6.5-8.5, as the addition of Link-A-Light Modifier provides the conditions necessary for efficient conjugation. Common non-buffering salts (e.g. sodium chloride), chelating agents (e.g. EDTA) and sugars may be present, as they have no effect on conjugation efficiency. Azide (0.02-0.1%) has little or no effect. Glycerol up to 50% has no effect. If the amine-free buffer is relatively concentrated and outside the pH range 6.6-8.5 you may need to add more Link-A-Light Modifier for each 10 µl of antibody. Excess Link-A-Light Modifier is provided so that you can check the pH of the buffer after the addition of the Modifier. Ideally the pH should be around 7.3-7.6, though efficient conjugation occurs anywhere between pH 6.8 and 7.8. Avoid buffer components that are nucleophilic, as these may react with Link-A-Light chemicals. Primary amines (e.g. amino acids, ethanolamine) and thiols (e.g. mercaptoethanol, DTT) fall within this class. Note: Unusually for an amine, Tris has little effect on conjugation efficiency as long as the concentration is 20mM or less.

Protocols:

1. **Before** you add antibody to the Link-A-Light Mix, add 1 µl of Link-A-Light Modifier reagent for each 10 µl of Antibody to be labeled. Mix gently.
2. Remove the screw cap from the vial of Link-A-Light Mix and pipette the antibody sample (with added Link-A-Light Modifier) directly onto the lyophilised material. Resuspend gently by withdrawing and re-dispensing the liquid once or twice using a pipette.
3. Place the cap back on the vial and leave the vial standing for 3 hours at room temperature (20-25°C). Alternatively, and sometimes more conveniently, conjugations can be set up and left overnight, as the longer incubation time has no negative effect on the conjugate.
4. After incubating for 3 hours (or more), add 1 µl of Link-A-Light Quencher FD reagent for every 10 µl of antibody used. The conjugate can be used after 30 minutes.

Storage:


Storage of Conjugates: For any new conjugate, initial storage at 4°C is recommended. A preservative may be desirable for Long-term storage. Other storage conditions (e.g. Frozen at -70°C or stored at -20°C with 50% glycerol) may also be satisfactory. The best conditions for any particular conjugate must be determined by experimentation.

Storage of Kit: The kit is shipped at ambient temperature in a tamper-evident polypropylene container. Store at -20°C upon receipt.

Pictures:

Absorbance Max (nm)	Emission Max (nm)	Extinction Coefficient (cm ⁻¹ M ⁻¹)	Fluorescent Colour	Stokes Shift
535	575	2 x 10 ⁶	Orange	40

Excitation (Left) and Emission (Right) Scan of Link-A-Light PE.

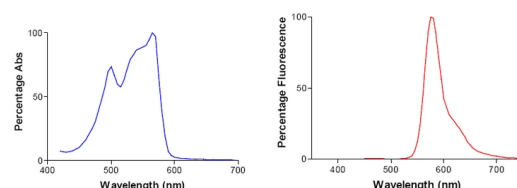

License:

Acris Antibodies' Link-A-Light conjugation kits are offered for research purposes alone, and are not intended for human, therapeutic or diagnostic use. The purchase of this conjugation kit conveys to the buyer (whether the buyer is a not-for-profit, academic or for-profit entity) the non-transferable right to use the amount of product purchased and the components of the product for in-house research. The buyer shall not sell or otherwise transfer this product, its components, or materials prepared therefrom to any third party. The buyer shall not use this product or its components for commercial purposes. For the avoidance of doubt, 'commercial purposes' means any activity by a party for consideration and includes, without limitation, use of the product or its components (i) in the manufacturing of conjugated materials (e.g. labeled antibodies), (ii) to provide a service, information or data, (iii) for therapeutic, diagnostic or prophylactic purposes, or (iv) for repackaging/resale, whether or not such product or its components are resold for use in research. The use of this product by the buyer constitutes agreement with the terms of this limited use label license for Link-A-Light products. For further information please contact Acris Antibodies GmbH, Schillerstraße 5, D-32052 Herford, Germany. Tel: +49-5221-34606-0, Fax +49-5221-34606-11, e-mail: info@acris-online.de