

Tel: 888.999.1371 or 858.642.2058

Fax: 858.642.2046 Web: www.ebioscience.com

Web: www.ebioscience.com
E-mail: info@ebioscience.com

Product Information

Contents: Fluorescein isothiocyanate (FITC) anti-mouse GM-

CSF

Catalog Number: 11-7331 Sizes: 25 tests, 50 ug, 100 ug

Formulation: Phosphate buffer pH 7.2,


150 mM NaCl, 0.09% NaN₃

Storage Conditions: Store at 4°C

DO NOT FREEZE.

LIGHT-SENSITIVE MATERIAL

Clone: MP1-22E9 Isotype: Rat IgG2a, κ

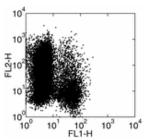

Mouse splenocytes were stimulated with ConA for 2 days, followed by IL-2 and IL-4 for 3 days, and restimulated with immobilized anti-CD3 and soluble anti-CD28 in the presence of Brefeldin A for 5 hours. The cells were surface stained with FITC anti-mouse CD4 (GK1.5) and intracellularly stained with PE anti-mouse GM-CSF (MP1-22E9).

Available	Formats of This Product			
Cat. No.	Format	Excite (nm)		Reported Applications
11-7331	FITC anti-mouse GM-CSF (Granulocyte/Macrophage-Colony Stimulating Factor, GMCSF)	488	518	IC
12-7331	PE anti-mouse GM-CSF (Granulocyte/Macrophage-Colony Stimulating Factor, GMCSF)	488	575	IC
14-7331	Affinity Purified anti-mouse GM-CSF (Granulocyte/Macrophage-Colony Stimulating Factor, GMCSF)	N/A	N/A	ELISA cap ELISPOT cap
16-7331	Functional Grade* Purified anti-mouse GM-CSF (Granulocyte/Macrophage-Colony Stimulating Factor, GMCSF)	N/A	N/A	ВА
*Function	al Grade™ (FG™): Azide-free, sterile-filtered, and endotoxin < 0.001 ng/μg.			

Functional Grade^{\mathbb{M}} (FG $^{\mathbb{M}}$): Azide-free, sterile-filtered, and endotoxin < 0.001 ng/ \mathbb{M} g. Purified: Contains azide, not sterile-filtered, and not endotoxin tested.

Description

The MP1-22E9 antibody reacts with mouse granulocyte/macrophage - colony stimulating factor (GM-CSF), a 14 kDa cytokine secreted by activated T cells and macrophages.

Usage

For research use only, not for diagnostic or therapeutic use. MP1-22E9 has been reported for use in intracellular flow cytometric analysis.

Applications Tested

The FITC conjugated MP1-22E9 is offered in 2 formats:

- μg size: has been tested by intracellular flow cytometric analysis and can be used at less than or equal to 0.5 μg/million cells. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.
- test size: has been pre-titrated and tested by intracellular flow cytometric analysis and can be used at 20 μl per million cells in a 100 μl total staining volume. The pre-titrated test size contains 0.2% BSA.

Related Products

 $\label{eq:copyright} \hbox{\mathbb{Q} 2000-2005 eBioscience, Inc.}$ Product For Research Use Only: Not for further distribution without written consent.