

OriGene Technologies, Inc.

9620 Medical Center Drive, Ste 200 Rockville, MD 20850 UNITED STATES Phone: +1-888-267-4436 Fax: +1-301-340-8606 techsupport@origene.com

OriGene Technologies GmbH

32052 Herford GERMANY Phone: +49-5221-34606-0 Fax: +49-5221-34606-11 info-de@origene.com

Schillerstr. 5

AM33124PU-T Monoclonal Antibody to basic Cytokeratin - Purified

Quantity: 20 μg Concentration: 0.2 mg/ml

Background: Cytokeratins are intermediate filament keratins found in the intracytoplasmic

cytoskeleton of epithelial tissue There are two types of Cytokeratins: the low weight, acidic type I cytokeratins and the high weight, basic or neutral type II. Cytokeratins are usually found in pairs comprising a type I Cytokeratin and a type II cytokeratin. The high molecular weight cytokeratins, which are the basic or neutral cytokeratins, comprise subtypes CK1 (67), CK2 (65.5), CK3 (64), CK4 (59), CK5 (58), CK6 (56), CK7 (54), CK8 (52.5) and CK9. The low molecular weight cytokeratins, which are the acidic cytokeratins, comprise subtypes CK10 (56.5), CK12 (56), CK13 (53), CK14 (50), CK16 (

48), CK17 (46), CK18 (45), CK19(48) and CK20(46).

Host / Isotype: Mouse / IgG1

Recommended Isotype

Controls:

SM10P (for use in human samples), SM20P (for use in rat samples), AM03095PU-N

Clone: 34betaE12

Immunogen: Solubilized Keratin extract from Human stratum corneum.

Format: State: Liquid purified IgG fraction from Bioreactor Concentrate

Purification: Protein A/G Chromatography

Buffer System: 10mM PBS

Preservatives: 0.05% Sodium Azide

Stabilizers: 0.05% BSA

Applications: Flow Cytometry: $0.5-1 \mu g/10^6$ cells.

Immunofluorescence: 0.5-1 µg/ml.

Western Blot: 0.5-1 µg/ml.

Immunohistochemistry on Frozen and Formalin-Fixed Paraffin Sections: $0.5-1\,\mu g/ml$

for 30 minutes at RT.

Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Citrate

Buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes.

Positive Control: PC12 cells, Skin, Prostate carcinoma.

Other applications not tested. Optimal dilutions are dependent on conditions and

should be determined by the user.

Molecular Weight: 67kDa (CK1); 58kDa (CK5); 56.5kDa (CK10); 50kDa (CK14)

Specificity:

This antibody recognizes CK1, CK5, CK10 and CK14. In normal epithelia, it stains stratified epithelia, myoepithelial cells and basal cells in the prostate gland and bronchi.

This Monoclonal Antibody shows no reactivity with hepatocytes, pancreatic acinar cells, proximal renal tubules, or endometrial glands; there is no reactivity with cells derived from simple epithelia. Mesenchymal tumors, lymphomas, melanomas, neural tumors, and neuroendocrine tumors are negative with this antibody. It stains myoepithelial cells and has been shown to be useful in distinguishing prostate adenocarcinoma from benign prostate.

This antibody has also been useful in separating benign from malignant intraductal breast proliferations.

Cellular Localization: Cytoplasmic. **Species:** Human, Mouse, Rat. Other species not tested.

Storage:

Store undiluted at 2-8°C.

Shelf life: one year from despatch.

General Readings:


- 1. Moinfar F, Man YG, Lininger RA, Bodian C, Tavassoli FA. Use of keratin 35betaE12 as an adjunct in the diagnosis of mammary intraepithelial neoplasia-ductal type--benign and malignant intraductal proliferations. Am J Surg Pathol. 1999 Sep;23(9):1048-58. PubMed PMID: 10478664.
- 2. Varma M, Linden MD, Amin MB. Effect of formalin fixation and epitope retrieval techniques on antibody 34betaE12 immunostaining of prostatic tissues. Mod Pathol. 1999 May;12(5):472-8. PubMed PMID: 10349984.

Pictures:

Formalin-Fixed, Paraffin-Embedded Human prostate (10X) stained

with Cytokeratin, HMW Antibody Cat.-No

AM33124PU (Clone 34bE12).

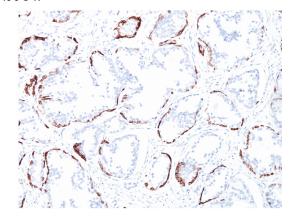



Formalin-Fixed, Paraffin-Embedded Human prostate (20X) stained with Cytokeratin, HMW Antibody Cat.-No AM33124PU (Clone 34bE12).

