

## OriGene Technologies, Inc.

9620 Medical Center Drive, Ste 200 Rockville, MD 20850 UNITED STATES Phone: +1-888-267-4436 Fax: +1-301-340-8606 techsupport@origene.com

## OriGene Technologies GmbH

Schillerstr. 5 32052 Herford GERMANY Phone: +49-5221-34606-0 Fax: +49-5221-34606-11 info-de@origene.com

## AM50121PU-T Monoclonal Antibody to Human Lambda Light Chain - Purified

**Quantity:** 20 μg

Concentration: 0.2 mg/ml

Host / Isotype: Mouse / IgG2a
Recommended Isotype AM03096PU-N

**Controls:** 

Clone: SPM559

Immunogen: Purified human IgG.

Genename: IGL

Format: State: Liquid purified IgG fraction from Bioreactor Concentrate

**Purification:** Protein A/G Chromatography

Buffer System: 10mM PBS

Preservatives: 0.05% Sodium Azide

Stabilizers: 0.05% BSA

**Applications: ELISA:** Use Antibody without BSA for Coating.

Flow Cytometry:  $0.5-1 \mu g/10^6$  cells. Immunofluorescence:  $1-2 \mu g/ml$ . Western Blot:  $0.5-1 \mu g/ml$ .

**Immunoprecipitation:** 0.5-1 μg/500 μg protein lysate.

Immunohistochemistry on Frozen and Formalin-Fixed Paraffin Sections: 0.5-1 µg/ml

for 30

minutes at RT.

Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Citrate

buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes.

*Positive Control*: 293T, Raji or hPBL cells. Tonsil or Spleen.

Other applications not tested. Optimal dilutions are dependent on conditions and

should be determined by the user.

Molecular Weight: ~22.5kDa

Specificity: This MAb is specific to lambda light chain of immunoglobulin and shows no cross-

reaction with lambda light chain or any of the five heavy chains. In mammals, the two light chains in an antibody are always identical, with only one type of light chain, kappa or lambda. The ratio of Kappa to Lambda is 70:30. However, with the

occurrence of multiple myeloma or other B-cell malignancies this ratio is disturbed. Antibody to the lambda light chain is reportedly useful in the identification of

leukemias, plasmacytomas, and certain non-Hodgkin's lymphomas. Demonstration of

clonality in lymphoid infiltrates indicates that the infiltrate is malignant.

*Cellular Localization*: Cell surface, cytoplasmic and secreted.

Species: Human.

Other species not tested.

Storage:

Store undiluted at 2-8°C.


Shelf life: one year from despatch.

**General Readings:** 

1. Campbell JP et. al. J Immunol Methods. 2013;391(1-2):1-13.

**Pictures:** 

Formalin-paraffin tonsil stained with Human Lambda Light Chain Antibody Cat.-No AM50121PU (Clone SPM559). Note cell membrane & cytoplasmic staining.

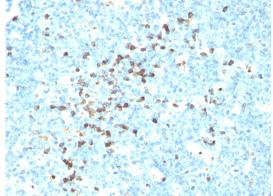

Western blot analysis of Lambda Light Chain in human intestine using Human Lambda Light Chain Antibody Cat.-No AM50121PU (Clone SPM559).

