

OriGene Technologies Inc.

9620 Medical Center Drive, Ste 200 Rockville, MD 20850 UNITED STATES Phone: +1-858-888-7900 Fax: +1-858-888-7904 <u>US-info@acris-antibodies.com</u>

BP7082 OriGene EU

Acris Antibodies GmbH Schillerstr. 5 32052 Herford GERMANY Phone: +49-5221-34606-0 Fax: +49-5221-34606-11 info@acris-antibodies.com

Polyclonal Antibody to 40S ribosomal protein S6 (RPS6) pSer236 - Aff - Purified BP7082 **Catalog No.: Quantity:** 0.1 ml **Background:** 40S ribosomal protein S6 (also known as RPS6) is a ~31 kDa substrate of p70 S6 kinase (p70S6K) and a major component of translational machinery involved in protein synthesis, cell growth, proliferation, and metabolism. Phosphorylation of RPS6 is rapamycin and wortmannin-sensitive as its activation is mediated by mTOR and PI3K pathways. Ribosomal protein S6 undergoes phosphorylation on multiple serines in the carboxyl terminal region in the order 236->235->240->244->247, due to the positions of these amino acid residues on the a-helix. Hyperphosphorylation of ribosomal protein S6 stimulates protein synthesis that mediates progression through the cell cycle. Host: Rabbit Chemically synthesized phosphopeptide derived from the region of human RPS6 that Immunogen: contains serine 236. **Remarks:** The sequence is conserved in mouse and rat. State: Liquid lg fraction Format: Purification: Sequential epitope-specific chromatography. The antibody has been negatively preadsorbed using a non-phosphopeptide corresponding to the site of phosphorylation to remove antibody that is reactive with non-phosphorylated RPS6. The final product is generated by affinity chromatography using a RPS6-derived peptide that is phosphorylated at serine 236. **Buffer System:** Dulbecco's phosphate buffered saline (without Mg2+ and Ca2+), pH 7.3 (+/- 0.1), 50% glycerol with 1.0 mg/mL BSA (lgG, protease free) as a carrier, containing 0.05 % sodium azide. Western blot (1:1000). **Applications:** Positive Control Used: HeLa +/- anisomycin. Other applications not tested. Optimal dilutions are dependent on conditions and should be determined by the user. **Specificity:** This antibody detects RPS6. It does not cross react with RPS6 phosphorylated on serines 244 and 247. Species: Human, Mouse, Rat. Other species not tested. Store the antibody at -20 °C. Can be shipped at 2 - 8 °C. Storage: Avoid repeated freezing and thawing. Centrifuge vial before opening. Shelf life: One year from despatch.

For research and in vitro use only. Not for diagnostic or therapeutic work. Material Safety Datasheets are available at www.acris-antibodies.com or on request.



OG/20121030

Acris Antibodies is now part of the OriGene family. Learn more at www.origene.com



BP7082: Polyclonal Antibody to 40S ribosomal protein S6 (RPS6) pSer236 - Aff -Purified

AN ORIGENE COMPANY	Purified
General Readings:	Pende, M., et al. (2004) S6K1(-/-)/S6K2(-/-) mice exhibit perinatal lethality and rapamycin- sensitive 5'-terminal oligopyrimidine mRNA translation and reveal a mitogen-activated protein kinase-dependent S6 kinase pathway. Mol. Cell. Biol. 24(8):3112-3124. Mourani, P.M., et al. (2004) Unique, highly proliferative growth phenotype expressed by embryonic and neointimal smooth muscle cells is driven by constitutive Akt, mTOR, and p70S6K signaling and is actively repressed by PTEN. Circulation 109(10):1299-1306. Lekmine, F., et al. (2004) Interferon-gamma engages the p70 S6 kinase to regulate phosphorylation of the 40S S6 ribosomal protein. Exp. Cell. Res. 295(1):173-182. Tuhackova, Z., et al. (2004) IL2-dependent phosphorylation of 40S ribosomal protein S6 is controlled by PI-3K/mTOR signalling in CTLL2 cells. Int. J. Mol. Med. 13(4):601-605. Shah, O.J., et al. (2003) Mitotic regulation of ribosomal S6 kinase 1 involves Ser/Thr, Pro phosphorylation of consensus and non-consensus sites by Cdc2. J. Biol. Chem. 278(18):16433-16442. Stewart, M.J. and G. Thomas (1994) Mitogenesis and protein synthesis: a role for ribosomal protein S6 phosphorylation? Bioessays 16(11):809-815. Review. Ferrari, S., et al. (1991) Mitogen-activated 70K S6 kinase. Identification of in vitro 40S ribosomal S6 phosphorylation sites. J. Biol. Chem. 266(33):22770-22775.
Protocols:	Western Blotting Procedure
	 Lyse approximately 10e7 cells in 0.5 mL of ice cold Cell Lysis Buffer (formulation provided below). This buffer, a modified RIPA buffer, is suitable for recovery of most proteins, including membrane receptors, cytoskeletal-associated proteins, and soluble proteins. Other cell lysis buffer formulations, such as Laemmli sample buffer and Triton-X 100 buffer, are also compatible with this procedure. Additional optimization of the cell stimulation protocol and cell lysis procedure may be required for each specific application. Remove the cellular debris by centrifuging the lysates at 14,000 x g for 10 minutes. Alternatively, lysates may be ultracentrifugedat 100,000 x g for 30 minutes for greater clarification. Carefully decant the clarified cell lysates into clean tubes and determine the protein concentration using a suitable method, such as the Bradford assay. Polypropylene tubes are recommended for storing cell lysates. React an aliquot of the lysate with an equal volume of 2x Laemmli Sample Buffer (125 mM Tris, pH 6.8, 10% glycerol, 10% SDS, 0.006% bromophenol blue, and 130 mM dithiothreitol [DTT]) and boil the mixture for 90 seconds at 100°C. Load 10-30 µg of the cell lysate into the wells of an appropriate single percentage or gradient minigel and resolve the proteins by SDS-PAGE. In preparation for the Western transfer, cut a piece of PVDF membrane slightly larger than the gel. Soak the membrane in methanol for 1 minute, then rinse with ddH2O for 5 minutes. Alternatively, nitrocellulose may be used. Soak the membrane into the sandwich apparatus. Transfer the proteins at 140 mA for 60-90 minutes at room temperature. Following the transfer, rinse the membrane with Tris buffered saline for 2 minutes. Block the membrane with blocking buffer (formulation provided below) for one hour at room temperature or overnight at 4°C. Lucubate the blocked blot with primary antibody at

13. Wash the blot with several changes of Tris buffered saline supplemented with 0.1% Tween 20.

14. Detect the antibody band using an appropriate secondary antibody, such as goat F(ab)2

For research and in vitro use only. Not for diagnostic or therapeutic work. Material Safety Datasheets are available at www.acris-antibodies.com or on request.

Acris Antibodies is now part of the OriGene family. Learn more at www.origene.com





anti-rabbit IgG alkaline phosphatase conjugate or goat F(ab)2 anti-rabbit IgG horseradish peroxidase conjugate in conjunction with your chemiluminescence reagents and instrumentation.

Cell Lysis Buffer Formulation: 10 mM Tris, pH 7.4 100 mM NaCl 1 mM EDTA 1 mM EGTA 1 mM NaF 20 mM Na4P207 2 mM Na3VO4 0.1% SDS 0.5% sodium deoxycholate 1% Triton-X 100 10% glycerol 1 mM PMSF (made from a 0.3 M stock in DMSO) or 1 mM AEBSF (water soluble version of PMSF) $60 \,\mu g/mL$ aprotinin $10 \mu g/mL$ leupeptin $1 \mu g/mL$ pepstatin (alternatively, protease inhibitor cocktail such as Sigma Cat. # P2714 may be used)

Transfer Buffer Formulation: 2.4 gm Tris base 14.2 gm glycine 200 mL methanol Q.S. to 1 liter, then add 1 mL 10% SDS. Cool to 4°C prior to use.

Tris Buffered Saline Formulation: 20 mM Tris-HCl, pH 7.4 0.9% NaCl

Blocking Buffer Formulation: 100 mL Tris buffered saline 5 gm BSA 0.1 mL Tween 20

Peptide Competition Experiment

The specificity of a Phosphorylation Site Specific Antibody (PSSA) in each experimental system can be confirmed through peptide competition. In this technique, aliquots of antibody are pre-incubated with peptide containing the sequence of the phosphopeptide immunogen used to raise the PSSA and the corresponding non-phosphopeptide. Following preincubation with the peptide, each antibody preparation is then used as a probe in antibody-based detection methods, such as Western blotting, immunocytochemistry, flow cytometry, or ELISA. With a PSSA specific for the phosphorylated target protein, pre-incubation with an excess of peptide containing the sequence of the phosphopeptide immunogen will block all antigen binding sites, while pre-incubation with the corresponding non-phosphopeptide will not affect the antibody.

For research and in vitro use only. Not for diagnostic or therapeutic work. Material Safety Datasheets are available at www.acris-antibodies.com or on request.





OG/20121030



In performing the Peptide Competition Experiment, it is important to note that the optimal dilutions of both antibody and peptide should be determined empirically for each specific application. The optimal dilution of antibody in these procedures is below saturating, as determined by previous experiments in your system.

The optimal dilution of peptide used in these procedures will depend on the overall affinity or avidity of the antibody, as well as the quantity of the target antigen. A 50-150 fold molar excess of peptide to antibody is found to be effective for most peptide competition experiments.

In the example presented below, the PSSA is used as a dilution of 1:1000 and the peptides are used at a concentration of 333 nM. The total volume of the phosphopeptide and nonphosphopeptide pre-incubated antibody preparations is 2 mL, sufficient for probing Western blot strips, as well as for use in other antibody-based detection methods. Under these conditions, the molar excess of peptide to antibody is > / = 50.

Procedure:

1. Prepare three identical test samples, such as identical PVDF or nitrocellulose strips to which the protein of interest has been transferred. The test samples should be blocked using a blocking buffer, such as Tris buffered saline supplemented with 0.1% Tween 20, and either 5% BSA or 5% non-fat dried milk.

2. Prepare 6.5 mL of working antibody stock solution (1:1000 in this example) by adding 6.5 µL of antibody stock solution to 6.5 mL of buffer containing blocking protein, such as TBS supplemented with 0.1% Tween 20, and either 3% BSA or 3% non-fat dried milk.

3. Apportion the unused PSSA into working aliquots and store at -20°C for future use (the stock PSSA contains 50% glycerol and will not freeze at this temperature).

4. Allow the lyophilized control peptides to reach room temperature, ideally under desiccation. 5. Reconstitute each of the control peptides to a concentration of 66.7 µM with nanopure water. (i.e. for a peptide with a molecular mass of 1500, reconstitution with 1 mL water yields a solution with a concentration of 66.7 μ M).

6. Apportion the unused reconstituted peptide solutions into working aliquots and store at -20°C for future use.

7. Label 3 test tubes as follows:

- tube 1: water only no peptide control

- tube 2: phosphopeptide

- tube 3: non-phosphopeptide

8. Into each tube, pipette the following components

- tube 1: 2 mL diluted PSSA solution plus 10 µL nanopure water

- tube 2: 2 mL diluted PSSA solution plus 10 µL phosphopeptide

- tube 3: 2 mL diluted PSSA solution plus 10 µL non-phosphopeptide

9. Incubate the three tubes for 30 minutes at room temperature with gentle rocking. During this incubation, the peptides have the chance to bind to the combining site of the antibody.

10. At the end of the incubation step, transfer the contents of each of the three tubes to clean reaction vessels containing one of the three identical test samples. For Western blotting strips:


Incubate the strips with the pre-incubated antibody preparations for 1 hour at room temperature or overnight at 4°C.

Wash each strip four times, five minutes each, to remove unbound antibody. Transfer each strip to a new solution containing a labeled secondary antibody [e.g., goat F(ab)2 antirabbit IgG alkaline phosphatase conjugate or goat F(ab)2 anti-rabbit IgG horseradish peroxidase conjugate.

Remove unbound secondary antibody by thorough washing, and develop the signal using your chemiluminescent reagents and instrumentation.

The signal obtained with antibody incubated with the "Water Only, No Peptide Control"

For research and in vitro use only. Not for diagnostic or therapeutic work.

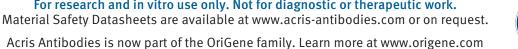



OG/20121030


4/5



(Tube 1), represents the maximum signal in the assay. This signal should be eliminated by preincubation with the "Phosphopeptide" (Tube 2), while pre-incubation with the "Non-Phosphopeptide" (Tube 3) should not impact the signal. If the "Phosphopeptide" only partially eliminates the signal, repeat the procedure using twice the volume of water or peptide solutions listed in Step 8. If partial competition is seen following pre-incubation with the "Non-Phosphopeptide", repeat the procedure using half the volumes of water or peptide solutions listed in Step 8.

Pictures:

Up-regulation and Antibody-Peptide **Competition** Extracts of HeLa cells unstimulated (1) or stimulated with 10 μ g/mL anisomycin for 10 minutes (2-6) were resolved by SDS-PAGE on a 10% Trisglycine gel and transferred to PVDF. The membrane was blocked with a 5% BSA-TBST buffer for one hour at room temperature, and incubated with ribosomal protein S6 [pS236] antibody for two hours at room temperature in a 3% BSA-TBST buffer, following prior incubation with: no peptide (1, 2), the phosphopeptide corresponding to ribosomal protein S6 [pSpS244/247] (3), the non-phosphopeptide corresponding to the phosphopeptide immunogen (4), a generic phosphoserine-containing peptide (5), or the phosphopeptide immunogen (6). After washing the membrane was incubated with goat F(ab')2 anti-rabbit IgG HRP conjugate and signals were detected using the Pierce SuperSignalTM method. The data show that only the phosphopeptide corresponding to ribosomal protein S6 [pS236] completely blocks the signal, demonstrating the specificity of the antibody. The data also show upregulation of the signal upon stimulation with anisomycin in this cell system.

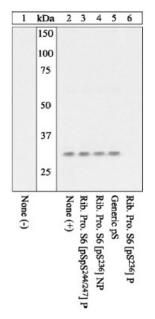

For research and in vitro use only. Not for diagnostic or therapeutic work. Material Safety Datasheets are available at www.acris-antibodies.com or on request. Acris Antibodies is now part of the OriGene family. Learn more at www.origene.com



OG/20121030