

## OriGene Technologies Inc.

9620 Medical Center Drive, Ste 200 Rockville, MD 20850 UNITED STATES

Phone: +1-888-267-4436 Fax: +1-301-340-8606 techsupport@origene.com

### OriGene EU

CI 077R

#### **Acris Antibodies GmbH**

Schillerstr. 5 32052 Herford GERMANY

Phone: +49-5221-34606-0 Fax: +49-5221-34606-11 info@acris-antibodies.com

# Monoclonal Antibody to T Cell Receptor (TCR) Vb 4 - PE

Alternate names: TCR V beta-4, TCR Vb4

Catalog No.:CL077RQuantity:50 μgConcentration:0.1 mg/mlHost / Isotype:Rat / IgG2aClone:CTVB4

Format: State: Liquid purified

Buffer System: PBS, 0.09% NaN3 and EIA grade BSA as a stabilizing protein to bring total

protein concentration to 4-5 mg/ml.

Label: PE

**Applications:** Flow cytometry.

Other applications not tested. Optimal dilutions are dependent on conditions and should

be determined by the user.

Specificity: This antigen monoclonal antibody reacts with theTCR Vb4 bearing T cells. TCR Vb4 is

expressed in most know mouse strains.

Species: Mouse.

Other species not tested.

**Storage:** Store the antibody undiluted at 2-8°C.

DO NOT FREEZE!

This product is photosensitive and should protected from light.

Shelf life: one year from despatch.

General Readings: 1. Tomonari, K., lovering, E. and Spencer, S. 1990 Correlation between the Vb4+ CD8+ T cell

population and the H-2d haplotype. Immunogenetics 31:333-339.

2. Padula, S.J., Lingenheld, E.G., Stabach, P.R., Chou, C.J., Kono, D.H. and Clark, R.B. 1991 Identification of encephalitogenic Vb4-bearing T cells in SJL mice-further evidence for the V

region disease hyopothesis? J. Immunol. 146:879-883.

## **Protocols:** FLOW CYTOMETRY ANALYSIS:

#### Method:

1. Prepare a cell suspension in media A. For cell preparations, deplete the red blood cell population with Lympholyte®-M cell separation medium.

2. Wash 2 times.

3. Resuspend the cells to a concentration of 2x10e7 cells/ml in media A. Add 50  $\mu$ l of this suspension to each tube (each tube will then contain 1 x 10e6 cells, representing 1 test).

4. To each tube, add ~1.0 μg\* of this Ab per 10e6 cells.

5. Vortex the tubes to ensure thorough mixing of antibody and cells.

TUV HORD CERT CHIEF TO SUIT OF THE PROPERTY CHIEF TO SUIT OF THE P



#### CL077R: Monoclonal Antibody to T Cell Receptor (TCR) Vb 4 - PE

- 6. Incubate the tubes for 30 minutes at 4°C. (It is recommended that the tubes are protected from light, since most flurochromes are light sensitive.)
- 7. Wash 2 times at 4°C.
- 8. Resuspend the cell pellet in 50 μl ice cold media B.
- 9. Transfer to suitable tubes for flow cytometric analysis containing 15  $\mu$ l of propidium iodide at 0.5 mg/ml in PBS. This stains dead cells by intercalating in DNA.

#### Media:

A. Phosphate buffered saline (pH 7.2) + 5% normal serum of host species + sodium azide (100  $\mu$ l of 2M sodium azide in 100 mls).

B. Phosphate buffered saline (pH 7.2) + 0.5% Bovine serum albumin + sodium azide (100  $\mu$ l of 2M sodium azide in 100 mls).

#### Results - Tissue Distribution by Flow Cytometry Analysis:


Mouse Strain: C3H.SW

Cell Concentration: 1x10e6 cells per tests
Antibody Concentration Used: 1.0 µg/10e6 cells

Isotypic Control: PE Rat IgG2a

**Pictures:** 

ΡΕ Vβ4

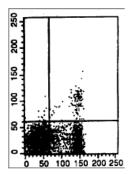

FITC CD4/CD8

Cell source: Spleen

Percentage of cells stained above control: 2.98%