

# NEUROMICS




## mGluR2/3

## Data Sheet

<b>Catalog Number:</b>	RA13102	<b>Host:</b>	Rabbit
<b>Product Type:</b>	Affinity Purified Antiserum	<b>Species Reactivity:</b>	Rat (only tested in rat)
<b>Immunogen Sequence:</b>	NGREVDSTSSSL Corresponding to the carboxy-terminus of rat mGluR2	<b>Format:</b>	Affinity purified. Sent in liquid form in Tris Glycine with 0.1% sodium azide and 1% BSA, pH 7.8.
<b>Applications:</b>	Immunohistochemistry 1:250 Immunocytochemistry 1:250 Western Blotting 1:1000 Dilutions listed only as a recommendation. Optimal dilution should be determined by investigator.		
<b>Storage:</b>	Store frozen. Aliquot as undiluted serum and immediately place at -20°C. Serum may have become trapped in top of vial during shipping. Centrifugation of vial is recommended before opening. Stable for at least 6 months at -20°C. Repeated freeze/thaw cycles compromise the integrity of the antiserum.		
<b>References:</b>	Tanabe, Y., Masu, M., Ishii, T., Shigemoto, R., and Nakanishi, S. (1992). <i>A family of metabotropic glutamate receptors</i> . <i>Neuron</i> 8 (1), 169-79.		

### Application Notes

Anti-mGluR2/3 affinity purified antisera can be used with standard protocols for immunohistochemistry, immunocytochemistry and Western Blotting. This antibody recognizes a conserved sequence in both metabotropic glutamate receptor 2 and metabotropic glutamate receptor 3.

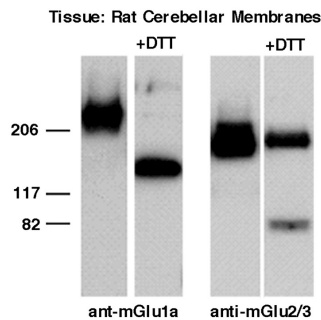

**Western blotting:** Western blot analysis of rat brain shows bands migrating at  $M_r=100,000$  and  $190,000$  which may be a dimer of the smaller band.

**Immunohistochemistry:** Antiserum was used on perfusion fixed tissue. Perfusion: 1) calcium-free Tyrode's solution, 2) paraformaldehyde-picric acid fixative, and 3) 10% sucrose in PBS as a cryo-protectant. Desired tissues were dissected and stored overnight in 10% sucrose in PBS.

Slide-mounted tissue sections were processed for indirect immunofluorescence. Slides were incubated with blocking buffer for 1 hour at room temperature. Primary antiserum was diluted with blocking buffer to the appropriate working concentration. Blocking buffer was removed and slides were incubated for 18-24 hours at 4°C with primary antiserum. Slides were rinsed 3 times and then incubated with secondary antibodies for 1 hour at room temperature. Slides were again rinsed 3 times and coverslipped. Staining was examined using fluorescence microscopy.

*Note:* Sodium azide ( $\text{NaN}_3$ ) interferes with peroxidase reactions and should not be used with peroxidase methodologies. If sodium azide is present in any steps of the staining procedure, the tissue should thoroughly be rinsed with sodium azide-free buffer before performing the peroxidase reaction.

### FOR RESEARCH USE ONLY

NEUROMICS REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

Please visit our website for complete information on protocols, buffers and solutions @ [www.neuromics.com](http://www.neuromics.com)

Neuromics, Inc. • 1313 5th Street S.E. Suite 113 • Minneapolis, MN 55414  
phone 877.327.3001 • phone 612.379.3604 • fax 612.379.3605